

Metadata of the chapter that will be visualized online

Chapter Title	Ratiometric Fluorescent Safranin-O Staining Allows the Quantification of Lignin Contents In Muro
Copyright Year	2022
Copyright Holder	The Author(s), under exclusive license to Springer Science+Business Media, LLC, part of Springer Nature
Author	Family Name Morel Particle Given Name Oriane Suffix Organization Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle Address Lille, France Division Institute of Biophysics Organization University of Natural Resources and Life Sciences Vienna Address Vienna, Austria
Author	Family Name Spriet Particle Given Name Corentin Suffix Organization Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle Address Lille, France Division Université de Lille, CNRS, Inserm, CHU Lille Organization Institut Pasteur de Lille, US 41-UMS 2014-PLBS Address Lille, France
Author	Family Name Lion Particle Given Name Cédric Suffix Organization Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle Address Lille, France
Author	Family Name Baldacci-Cresp Particle Given Name Fabien Suffix Organization Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle

	Address	Lille, France
Author	Family Name	Pontier
	Particle	
	Given Name	Garance
	Suffix	
	Organization	Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle
	Address	Lille, France
Author	Family Name	Baucher
	Particle	
	Given Name	Marie
	Suffix	
	Division	Laboratoire de Biotechnologie Végétale (LBV)
	Organization	Université Libre de Bruxelles
	Address	Gosselies, Belgium
Author	Family Name	Biot
	Particle	
	Given Name	Christophe
	Suffix	
	Organization	Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle
	Address	Lille, France
Author	Family Name	Hawkins
	Particle	
	Given Name	Simon
	Suffix	
	Organization	Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle
	Address	Lille, France
Corresponding Author	Family Name	Neutelings
	Particle	
	Given Name	Godfrey
	Suffix	
	Organization	Université de Lille, CNRS, UMR 8576, UGSF – Unité de Glycobiologie Structurale et Fonctionnelle
	Address	Lille, France
	Email	godfrey.neutelings@univ-lille.fr

Abstract In some specific vascular plant tissues, lignin can impregnate the entire cell wall to make it more rigid and hydrophobic. Different techniques have been developed in the past years to make possible the quantification of this

polyphenolic polymer at the organ or tissue level, but difficulties of access to the cellular level remain. Here we describe an approach based on ratiometric emission measurements using safranin-O and the development of a macro adapted for the FIJI software, which makes it possible to quantify lignin in three different layers of the cell wall on images captured on a fluorescent confocal microscope.

Keywords (separated by '-')	Lignin quantification - Safranin-O - Cell wall - Segmentation - Confocal microscopy
--------------------------------	---

Ratiometric Fluorescent Safranin-O Staining Allows the Quantification of Lignin Contents In Muro

Oriane Morel, Corentin Spriet, Cédric Lion, Fabien Baldacci-Cresp, Garance Pontier, Marie Baucher, Christophe Biot, Simon Hawkins, and Godfrey Neutelings

AU1

Abstract

In some specific vascular plant tissues, lignin can impregnate the entire cell wall to make it more rigid and hydrophobic. Different techniques have been developed in the past years to make possible the quantification of this polyphenolic polymer at the organ or tissue level, but difficulties of access to the cellular level remain. Here we describe an approach based on ratiometric emission measurements using safranin-O and the development of a macro adapted for the FIJI software, which makes it possible to quantify lignin in three different layers of the cell wall on images captured on a fluorescent confocal microscope.

Key words Lignin quantification, Safranin-O, Cell wall, Segmentation, Confocal microscopy

1 Introduction

Lignin is a major component of plant cell walls present in some specialized tissues of vascular plants. This polyphenolic polymer is mainly derived from the polymerization of *p*-coumaryl, coniferyl, and sinapyl alcohols differing in their degree of methoxylation. After their incorporation, these three *p*-hydroxycinnamyl alcohol precursors (monolignols) give rise to *p*-hydroxyphenyl (H), guaiacyl (G), and syringyl (S) units, respectively. Lignin is responsible for water transport in the xylem vessel elements and confers strength to cell walls contributing to the general mechanical support necessary for plant upright growth [1]. Studies on lignin biosynthesis are motivated by the need to understand its role in plant physiology but also by the increasing use of biomass in agro-industry [2].

A good understanding of lignin metabolism necessarily requires the use of techniques allowing its quantification. Two main “wet chemistry” approaches are commonly used for

quantifying lignin, depending upon the plant species and, more importantly, whether the tissues contain more or less quantities of lignins. Thus, the method using the solubilization of lignins by acetyl bromide [3, 4] is well adapted to herbaceous species whereas the gravimetric Klason method [5] is more suitable for woody species. There are also different methods to determine the H, G, and S subunit composition and the S/G ratio of lignins. Historically, thioacidolysis [6] and nitrobenzene oxidation [7] methods were used, but they are now gradually replaced by pyrolysis coupled with gas chromatography/mass spectrometry (GC/MS) [8] or nuclear magnetic resonance (NMR) [9]. These approaches provide relevant information on the different tissues/organs studied, but their use leads to irreversible tissue destruction. It is therefore not possible to obtain information concerning lignin quantities and compositions at the cellular level.

Therefore, histochemical approaches were also developed such as phloroglucinol staining (Wiesner reaction), which specifically reveals the cinnamaldehyde functions of S and G unit derivatives. In this case, the intensity of the coloration will depend on the quantity of the red chromogen produced by the reaction [10]. The Maïle test, on the other hand, uses potassium permanganate and hydrochloric acid to transform G and S residues into catechols, which then give rise to orange-brown quinones in the case of G lignins and purple-red quinones for G-S lignins [11]. Other approaches exploit the autofluorescence of lignins [12] or the possibility to visualize them by using fluorescent compounds such as auramine [13] and acriflavine [14]. All these techniques can therefore provide relatively detailed spatial information, but without offering the possibility to produce quantitative data.

Safranin-O has an interesting potential as a fluorescent probe [15]. Indeed, the emission spectrum of safranin-O depends on the amount of lignin present in the cell wall (Fig. 1). When present in high amounts, the emission will be more intense in the red part. The emission spectrum (at an excitation of 488 nm) always shows a maximum at 568 nm. Before this value, the aspect of the graph is invariable whatever the amount of lignin. On the other hand, in the orange/red area between 568 and 600 nm, it is dependent on the cell wall environment. Indeed, emission values are higher in the lignified tissues (xylem) compared to the non-lignified tissues. Based on these observations, our group developed an approach including a ratiometric method allowing the relative quantification of lignins by comparing the emissions between 540–560 nm (stable portion) and 570–600 nm (variable portion) [16].

In this chapter, we first describe the methodology allowing the quantification of lignin at the cellular level following confocal microscopy acquisition. In a second step, we demonstrate the possibility to apply a segmentation approach [17] to subdivide the cell wall on these images in an automatic way in order to obtain

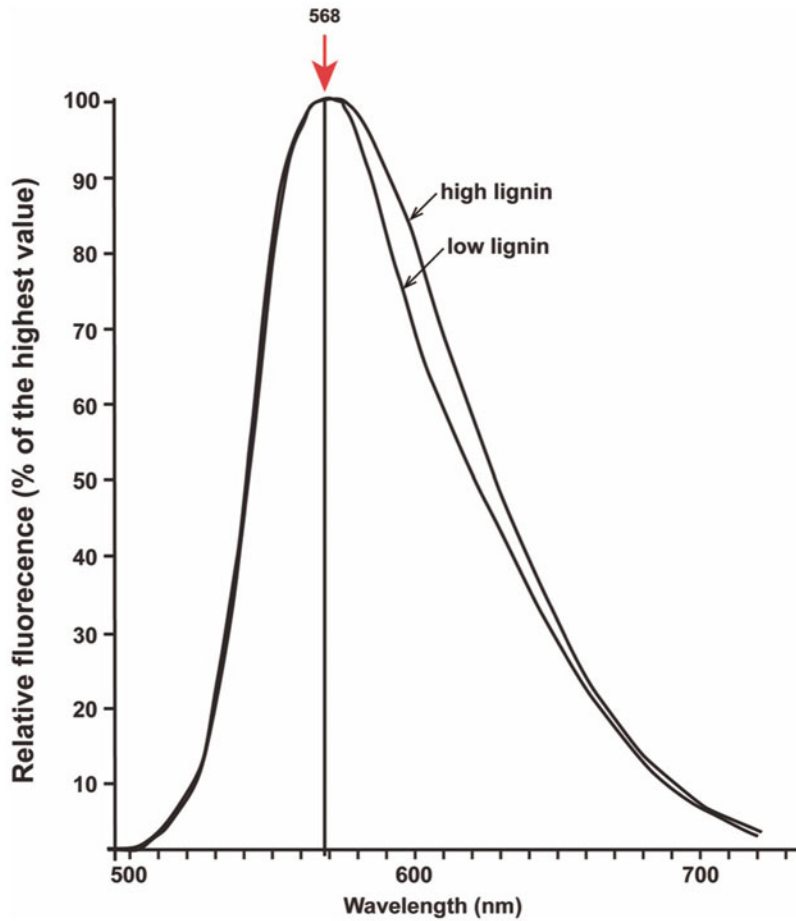


Fig. 1 Average safranin-O spectra (with excitation at 488 nm) of cell walls from different flax (*Linum usitatissimum*) tissues containing low amounts (bast fiber and cambial cell walls) or high amounts (xylem cell walls) of lignins

information on lignification in the cell corners, the compound middle lamella, and the secondary cell walls.

2 Materials

2.1 Plant Growth

1. Plant seeds. In this study we used *Arabidopsis thaliana* as an example.
2. Pots filled with potting soil.

2.2 Sample Preparation

1. 50% and 70% (v/v) ethanol solutions.
2. 4% (w/v) agarose solution in ultrapure water.
3. Cyanoacrylate glue.

	4. 0.2% safranin-O in a 50% (v/v) aqueous ethanol solution kept in the dark.	88 89
2.3 Equipment and Software	1. Plant growth chamber in which temperature and artificial light variations are possible.	90 91
	2. Vibratome (e.g., VT-1000S, Leica, Wetzlar, Germany) for cutting sections to a thickness of 80 μm (<i>see Note 1</i>).	92 93
	3. Microplate shaker.	94
	4. Confocal laser scanning microscope [we use a Nikon A1R confocal (Nikon, Tokyo, Japan) equipped with a 60 \times /1.4 aperture oil immersion objective (Plan APO VC) and the NIS Element AR3.0 software].	95 96 97 98
	5. FIJI v2.0 [18]. This software can be downloaded at https://fiji.sc .	99 100
	6. The associated FIJI macro developed for the quantification of fluorescence and segmentation has the following assigned doi, 10.5281/zenodo.5775647, and can be downloaded from the Zenodo website (https://zenodo.org/record/5775647#.YbiaNi17RhE).	101 102 103 104 105
<hr/>		
3 Methods		106
3.1 Plant Growth	In our experiment, we used <i>Arabidopsis thaliana</i> seeds that were stratified during 72 h at 4 $^{\circ}\text{C}$ in the dark. Plants were cultivated in a growth chamber under 12 h light (100 PAR light intensity) at 22 $^{\circ}\text{C}$ and 12 h dark at 20 $^{\circ}\text{C}$ for 6 weeks. The photoperiod was then changed to 16 h light and 8 h dark, and the floral stems were cut when they reached a size of 30 cm.	107 108 109 110 111 112
3.2 Sample Preparation	1. Collect floral stem samples of 0.5-cm length.	113
	2. Melt the 4% agarose solution and maintain it at 55 $^{\circ}\text{C}$ (<i>see Note 2</i>).	114 115
	3. Pour 1.5 mL of this solution in enough wells of a 48-well cell culture plate, and let the agarose cool down until it becomes thicker but still stays liquid.	116 117 118
	4. Place the floral stem samples in the center of each well with tweezers, and maintain in a vertical position until the agarose hardens.	119 120 121
	5. Leave the plate on the bench for another 10 min, and then unmold with a thin spatula.	122 123
	6. Cut the basis of the agarose perfectly perpendicular to the fragment with a razor blade.	124 125
	7. Glue the agarose-embedded fragment on the specimen holder with cyanoacrylate glue (<i>see Note 3</i>).	126 127

8. Cut 80- μm -thick fragments with the vibratome (*see Note 4*). 128
 9. Place the sections in a 48-well cell culture plate containing 1-mL 50% (v/v) ethanol solution. 129
130
- 3.3 Sample Staining**
1. Replace the ethanol in the section-containing well by a fresh or daily thawed 0.2% safranin-O solution (*see Note 5*), and incubate for 10 min on a microplate shaker in the dark, at room temperature. The shaking intensity should be high enough to allow good mobility of the section in the liquid. 131
132
133
134
135
 2. Remove all safranin-O solution, and replace with 50% ethanol. Incubate with shaking for 10 min in the dark. 136
137
 3. Remove the solution and rinse twice with ultrapure water for 15 min under the same shaking parameters. 138
139
 4. Mount the sections in water on a glass slide and place a coverslip (*see Note 6*). 140
141
- 3.4 Image Capture**
1. Perform the acquisition on the tissues by using three independent channels. The first is on autofluorescence excitation at 405 nm and emission at 450/500 nm; the second is on the green channel with excitation at 488 nm and emission at 530–560 nm. The third is on the red channel with excitation at 561 nm and emission at 570–600 nm. All tracks must be acquired sequentially, and parameters must be kept constant for all acquisitions. 142
143
144
145
146
147
148
149
 2. Select a flat region and adjust your focus to perform images between 6 and 9 μm deep from the surface (*see Note 7*). 150
151
- 3.5 Image Analysis**
1. Crop the sharp area of the images acquired in Subheading 3.4. 152
 2. Treat images using FIJI v2.0 and the associated homemade macro developed for this method (*see Note 8*). First open the Fiji app; open the macro and finally the image cropped in 1. Run the macro and follow the instructions. The macro includes two steps: 153
154
155
156
157
 - (a) First, a ratio image is generated as previously described [16] (Fig. 2). The result appears as a ratiometric image with a color range between 0 (purple) and 2 (red). The intensity ratio values also appear in the results window and can be used for quantification. 158
159
160
161
162
 - (b) Combining safranin fluorescence and lignin autofluorescence information, the ratio image obtained in 3.5.1.1 is then segmented, as previously described [17], into three cell walls regions: the cell corner (CC), the compound middle lamella (CML) consisting of the middle lamella and the primary cell wall, and the secondary cell wall (SCW) (Fig. 3) (*see Note 9*). 163
164 AU2
165
166
167
168
169

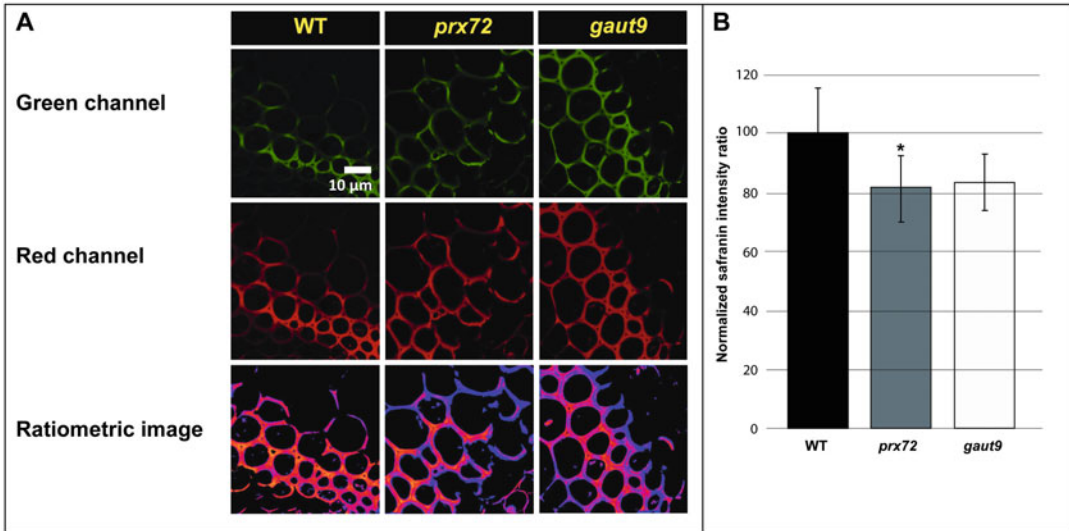


Fig. 2 Local lignin content determination on *Arabidopsis* floral stems (wild type and mutants) using safranin-0 fluorescence-based method. (a) Green and red channels correspond to safranin fluorescence. The ratiometric image was obtained by applying the specific FIJI macro on both images. (b) The normalized safranin intensity ratios (\pm SD) are represented (WT = 100%). The quantification was performed on interfascicular fibers. Means (\pm SD) with an asterisk are significantly different according to Student *t*-test (p -value < 0.05, $n = 3$). WT, wild type; *prx72*, Salk_136893 mutant; *gaut9*, SALK_040287 mutant

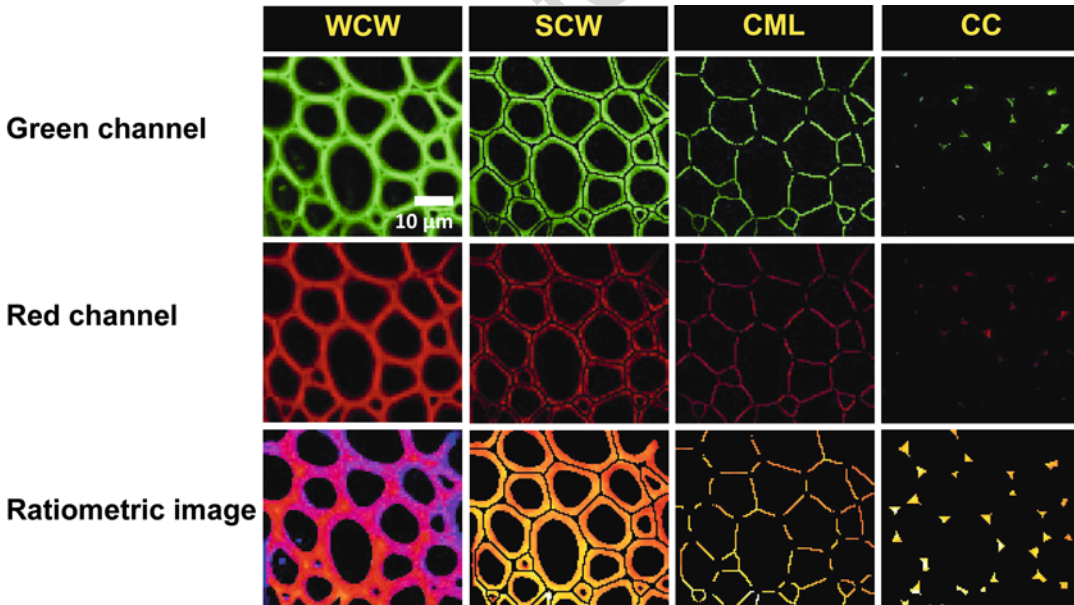


Fig. 3 Automatic segmentation of different cell wall zones. Three masks were applied to images of *Arabidopsis* flowering stem cross sections stained with safranin-0 and observed by confocal laser scanning microscopy. Green and red channels correspond to safranin-0 fluorescence. WCW, whole cell walls; SCW, secondary cell walls; CML, compound middle lamella; CC, cell corner

4 Notes

170

1. The use of a vibratome ensures regular homogenous sectioning of the plant tissues. However, hand sectioning with a sharp razor blade can also be performed. The sections should then be checked under the microscope, and only the most fine and regular ones should be selected for the safranin-O staining. 171-175
2. If using a vibratome, the agarose solution must be used within the week. In our protocol, we suggest to use multi-well cell culture plates, but it is also possible to use 1.5-mL microtubes. Their bottom conical part should first be cut off and placed inside the rest of the tube. Both parts are then sealed with a piece of parafilm. The melded agarose can then be poured in the tube and the sample placed in the center. When the gel has cooled down, it can be unmolded by pressing on the bottom conical portion of the tube. 176-184
3. On the Leica VT-1000S vibroslicer, it is possible to place up to three samples on one holder. 185-186
4. The speed and the frequency depend on the tissue. For a fragment taken at the base of the floral stem, use a speed of 0.4 mm/s (position 6 on the Leica model) and a frequency of 70 Hz (position 7). 187-190
5. The 0.2% safranin-O in a 50% (v/v) ethanol solution can be stored for several months at -20°C . When the solution is thawed, it should be stored in the dark and used within 24 h. 191-193
6. The samples must be stained and processed in the same day in order to be compared. Each characterization of mutants, for example, must be done with the corresponding wild type. 194-196
7. When using the confocal microscope, the focus is more accurate when using the red channel. 197-198
8. Accurate analysis should be performed on at least three different plants on which five sections were performed and at least two images captured per section (total = 30 images). 199-201
9. The segmentation macro also allows access to the intensity ratio values in the three layers of the cell wall. It is therefore possible to perform accurate quantification and statistics on lignin amounts. 202-205

207 References

- | | |
|---|--|
| <p>209 1. Boerjan W, Ralph J, Baucher M (2003) Lignin
210 biosynthesis. <i>Annu Rev Plant Biol</i> 54:
211 519–546. https://doi.org/10.1146/annurev.
212 arplant.54.031902.134938</p> | <p>2. Li X, Weng JK, Chapple C (2008) Improve- 213
ment of biomass through lignin modification. 214
<i>Plant J</i> 54:569–581. https://doi.org/10.
215 1111/j.1365-313X.2008.03457.x
216</p> |
|---|--|

- 217 3. Morrison I (1972) Improvements in the acetyl
218 bromide technique to determine lignin and
219 digestibility and its application to legumes. *J*
220 *Sci Food Agric* 23:1463–1469
- 221 4. Morrison I (1972) A semi-micro method for
222 the determination of lignin and its use in pre-
223 dicting the digestibility of forage crops. *J Sci*
224 *Food Agric* 23:455–463
- 225 5. Johnson D, Moore WE, Zank LC (1961) The
226 spectrophotometric determination of lignin in
227 small wood samples. *J Tech Assoc Pulp Pap Ind*
228 44:793–798
- 229 6. Lapierre C, Monties B, Rolando C (1986)
230 Thioacidolyse von Pappelholzligninen: Nach-
231 weis von monomeren Syringylverbindungen
232 und Charakterisierung von Guajacyl-Syringyl-
233 Ligninfraktionen. *Holzforschung* 40:
234 113–118. [https://doi.org/10.1515/hfsg.](https://doi.org/10.1515/hfsg.1986.40.2.113)
235 [1986.40.2.113](https://doi.org/10.1515/hfsg.1986.40.2.113)
- 236 7. Monties B (1989) Lignins. In: Dey PM, Har-
237 borne JB (eds) *Methods in plant biochemistry*.
238 Academic Press, pp 113–147
- 239 8. Wagner A, Ralph J, Akiyama T, Flint H,
240 Phillips L, Torr K et al (2007) Exploring ligni-
241 fication in conifers by silencing
242 hydroxycinnamoyl-CoA: shikimate hydroxy-
243 cinnamoyltransferase in *Pinus radiata*. *Proc*
244 *Natl Acad Sci U S A* 104:11856–11861.
245 <https://doi.org/10.1073/pnas.0701428104>
- 246 9. Mansfield SD, Kim H, Lu F, Ralph J (2012)
247 Whole plant cell wall characterization using
248 solution-state 2D NMR. *Nat Protoc* 7:
249 1579–1589. [https://doi.org/10.1038/nprot.](https://doi.org/10.1038/nprot.2012.064)
250 [2012.064](https://doi.org/10.1038/nprot.2012.064)
- 251 10. Clifford MN (1974) Specificity of acidic phlor-
252 oglucinol reagents. *J Chromatogr A* 94:
253 321–324. [https://doi.org/10.1016/S0021-](https://doi.org/10.1016/S0021-9673(01)92389-1)
254 [9673\(01\)92389-1](https://doi.org/10.1016/S0021-9673(01)92389-1)
- 255 11. Day A, Ruel K, Neutelings G, Crônier D,
256 David H, Hawkins S et al (2005) Lignification
257 in the flax stem: evidence for an unusual lignin
in bast fibers. *Planta* 222:234–245. [https://](https://doi.org/10.1007/s00425-005-1537-1) 258
doi.org/10.1007/s00425-005-1537-1 259
12. Donaldson LA, Radotic K (2013) Fluorescence 260
lifetime imaging of lignin autofluorescence in 261
normal and compression wood. *J Microsc* 251:
262 178–187. [https://doi.org/10.1111/jmi.](https://doi.org/10.1111/jmi.12059) 263
[12059](https://doi.org/10.1111/jmi.12059) 264
13. Pesquet E, Ranocha P, Legay S, Digonnet C, 265
Barbier O, Pichon M et al (2005) Novel mar- 266
kers of xylogenesis in zinnia are differentially 267
regulated by auxin and cytokinin. *Plant Physiol* 268
139:1821–1839. [https://doi.org/10.1104/](https://doi.org/10.1104/pp.105.064337) 269
[pp.105.064337](https://doi.org/10.1104/pp.105.064337) 270
14. Donaldson LA (2001) Lignification and lignin 271
topochemistry – an ultrastructural view. *Phyto-* 272
chemistry 57:859–873. [https://doi.org/10.](https://doi.org/10.1016/S0031-9422(01)00049-8) 273
[1016/S0031-9422\(01\)00049-8](https://doi.org/10.1016/S0031-9422(01)00049-8) 274
15. Bond J, Donaldson L, Hill S, Hitchcock K 275
(2008) Safranin fluorescent staining of wood 276
cell walls. *Biotech Histochem* 83:161–171. 277
[https://doi.org/10.1080/](https://doi.org/10.1080/10520290802373354) 278
[10520290802373354](https://doi.org/10.1080/10520290802373354) 279
16. Baldacci-Cresp F, Spriet C, Twyffels L, Bler- 280
vacq A-S, Neutelings G, Baucher M et al 281
(2020) A rapid and quantitative safranin- 282
based fluorescent microscopy method to evalu- 283
ate cell wall lignification. *Plant J* 102: 284
1074–1089. [https://doi.org/10.1111/tpj.](https://doi.org/10.1111/tpj.14675) 285
[14675](https://doi.org/10.1111/tpj.14675) 286
17. Morel O, Lion C, Neutelings G, Stefanov J, 287
Baldacci-Cresp F, Simon C et al (2022) 288
REPRISAL: mapping lignification dynamics 289
using chemistry, data segmentation, and ratio- 290
metric analysis. *Plant Physiol* 188(2):816–830. 291
<https://doi.org/10.1093/plphys/kiab490> 292
18. Schindelin J, Arganda-Carreras I, Frise E, 293
Kaynig V, Longair M, Pietzsch T et al 294
(2012) Fiji: an open-source platform for 295
biological-image analysis. *Nat Methods* 9: 296
676–682. [https://doi.org/10.1038/nmeth.](https://doi.org/10.1038/nmeth.2019) 297
[2019](https://doi.org/10.1038/nmeth.2019) 298

Author Queries

Chapter No.: 21 513975_2_En

Query Refs.	Details Required	Author's response
AU1	Please check whether the author names and affiliations are correct.	
AU2	Please provide appropriate cross-references instead of 3.5.1.1.	

Uncorrected Proof